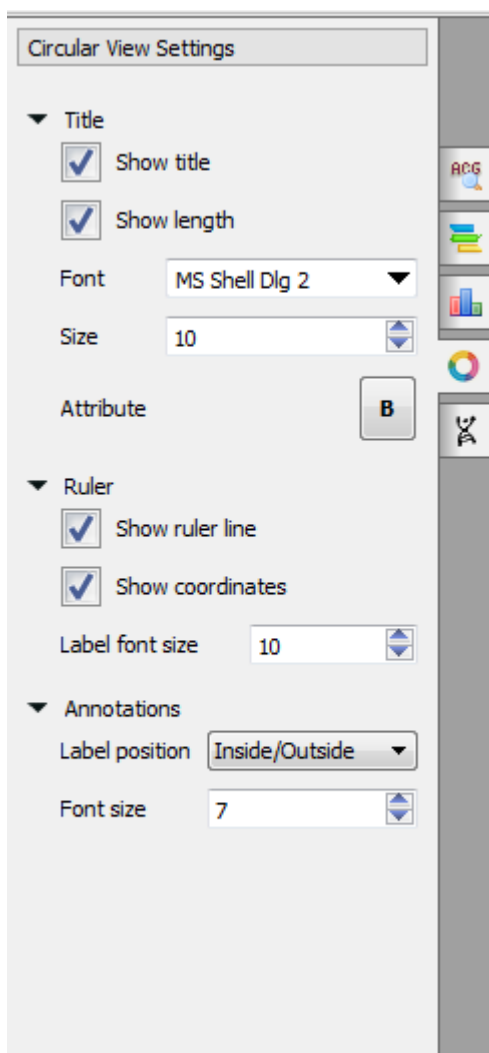


Circular View Settings

To configure circular view settings go to the *Circular View Settings* tab in the *Options Panel*.

Activate the circular view for a sequence and the following settings will appear:

The image shows a software dialog box titled "Circular View Settings". It is organized into three main sections, each with a downward-pointing triangle icon to its left. The first section, "Title", contains two checked checkboxes for "Show title" and "Show length", a "Font" dropdown menu set to "MS Shell Dlg 2", a "Size" spinner box set to "10", and an "Attribute" button with a bold "B" icon. The second section, "Ruler", contains two checked checkboxes for "Show ruler line" and "Show coordinates", and a "Label font size" spinner box set to "10". The third section, "Annotations", contains a "Label position" dropdown menu set to "Inside/Outside" and a "Font size" spinner box set to "7". To the right of the settings area is a vertical toolbar with icons for "ACG", a stack of papers, a bar chart, a color wheel, and a DNA double helix.

In the title section you can show or hide title and length, change font, size and attribute.

In the ruler section you can show or hide ruler line and coordinates and change the label font size.

In the annotation section you can select the label position and change the label size.

The following label positions are available:

- inside - all labels are inside of the annotations
- outside - all labels are outside of the annotations
- inside/outside - if the label can fit the annotation and it is not auto-annotation, it's located inside. Otherwise outside.
- none - no labels at all