

Tandem Repeats Finding

To find tandem repeats, select the *Analyze Find tandems...* context menu item in the *Sequence View* window.

In the opened dialog you can specify the tandem search parameters, the region to search in and the result parameters:

Find Tandems

Base | **Advanced**

Tandem finder parameters

Tandem preset: All

Min period: 1 n

Max period: 1000000 n

Region to process

Region: Whole sequence | 1 - 140426

Save annotation(s) to

☒ Existing table: NC_014267 features [NC_014267.1.gb]

☐ Create new table

☐ Use auto-annotations table

Annotation parameters

Group name: <auto>


Annotation name: repeat_unit

Description:

Start Cancel Help

The dialog parameters:

Tandem preset — specify the tandem repeats parameters with predefined values by selecting the available preset:

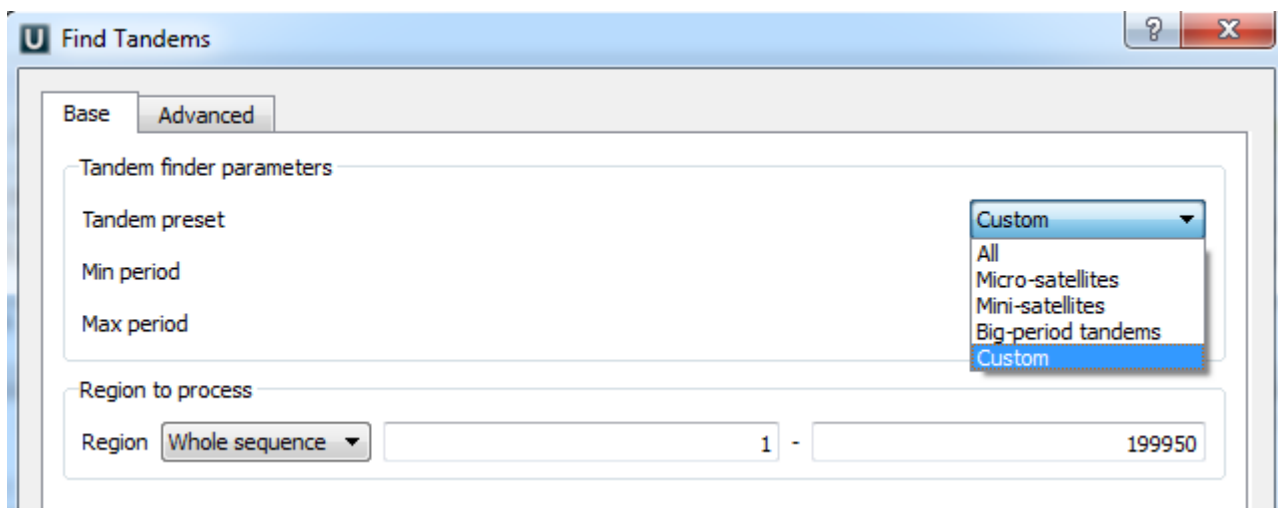


Min period, Max period — the minimum and maximum acceptable repeat length measured in base symbols.

Region to process — specify the region to search in the whole sequence, a custom region or the region of the current selection (if any).

In the *Save annotation(s) to group* you can set up a file to store annotations. It could be either an existing annotation table object, a new annotation table or auto-annotations table (if it is available).

In the *Annotation parameters* group you can specify the name of the group and the name of the annotation. If the group name is set to <auto> UGENE will use the group name as the name for the group. You can use the '/' characters in this field as a group name separator to create subgroups. If the annotation name is set to *by type* UGENE will use the annotation type from the *Annotation type*: table as the name for the annotation. Also you can add a description in the corresponding text field.

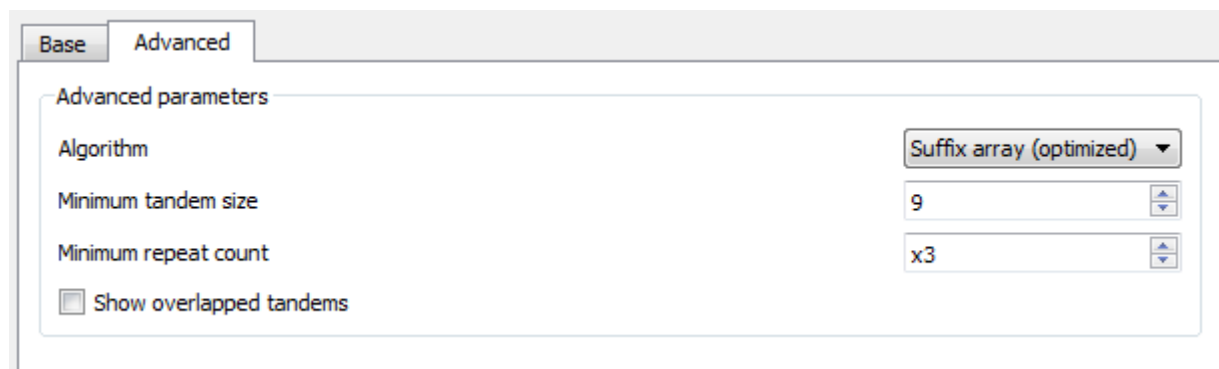

Additional search options can be found in the *Advanced* tab:

Algorithm — the algorithm parameter allows to select the search algorithm. The default and a fast one is optimized suffix array algorithm.

Minimum tandem size — the minimum tandem size sets the limit on minimum acceptable length of the tandem, i.e. the minimum total repeats length of the searched tandem.

Minimum repeat count — the minimum number of repeats of a searched tandem.

Show overlapped tandems — check if the plugin should search for the overlapped tandems, otherwise keep unchecked.

- Tandem Repeats Search Result