


Circular View Settings

To configure circular view settings go to the *Circular View Settings* tab in the *Options Panel*.

Activate the circular view for a sequence and the following settings will appear:

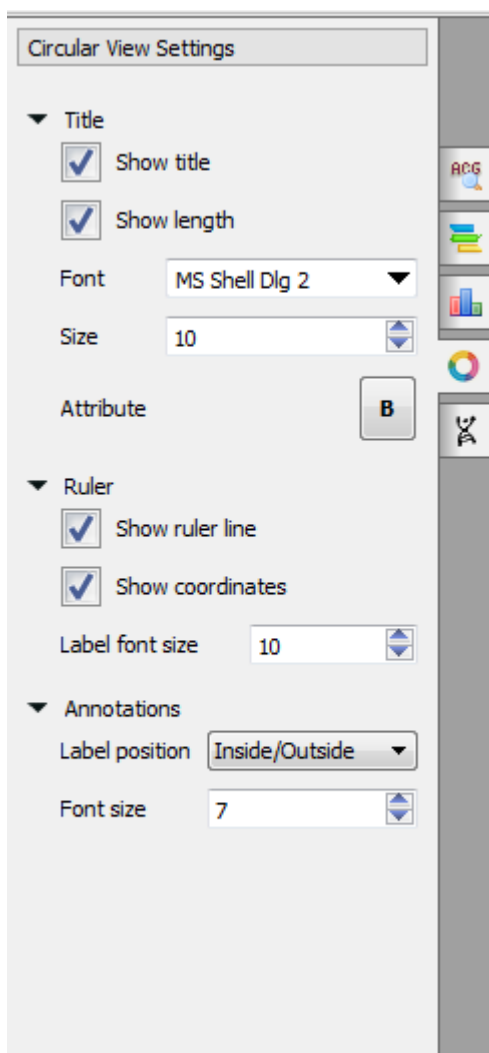

The image shows a software dialog box titled "Circular View Settings". It is organized into three main sections, each with a dropdown arrow to its left:

- Title**: Contains two checked checkboxes for "Show title" and "Show length". Below these are controls for "Font" (a dropdown menu showing "MS Shell Dlg 2") and "Size" (a numeric input field showing "10"). There is also an "Attribute" label next to a button with a bold "B" icon.
- Ruler**: Contains two checked checkboxes for "Show ruler line" and "Show coordinates". Below these is a "Label font size" control with a numeric input field showing "10".
- Annotations**: Contains a "Label position" dropdown menu showing "Inside/Outside" and a "Font size" control with a numeric input field showing "7".

On the right side of the dialog, there is a vertical toolbar with icons for "ACG", a stack of papers, a bar chart, a color wheel, and a DNA double helix.

In the title section you can show or hide title and length, change font, size and attribute.

In the ruler section you can show or hide ruler line and coordinates and change the label font size.

In the annotation section you can select the label position and change the label size.

The following label positions are available:

- inside - all labels are inside of the annotations
- outside - all labels are outside of the annotations
- inside/outside - if the label can fit the annotation and it is not auto-annotation, it's located inside. Otherwise outside.
- none - no labels at all