


Export Chromatogram Consensus

To export consensus sequence use the *Consensus* tab of the *Options Panel*:

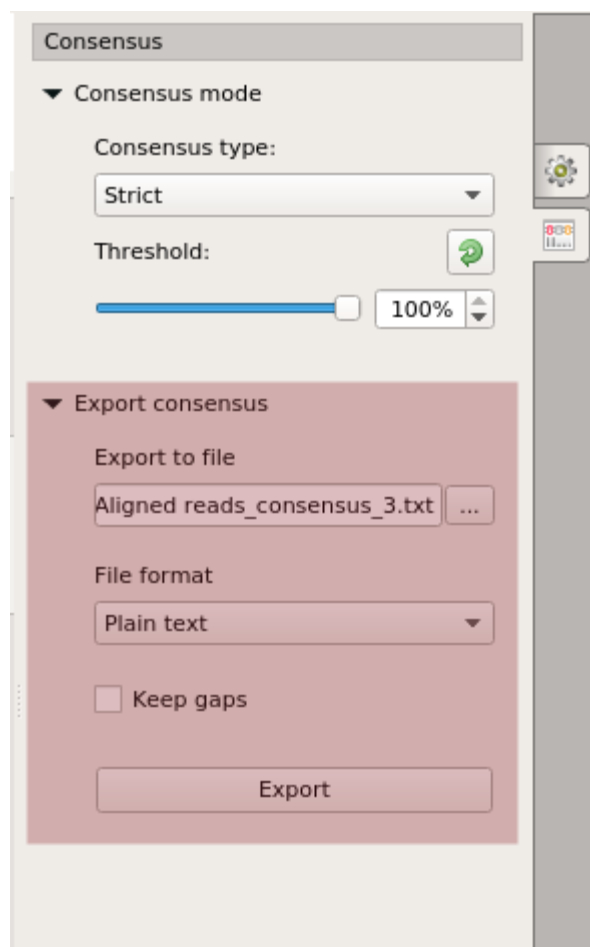

The image shows a software interface for exporting a consensus sequence. It features a 'Consensus' tab at the top. Below the tab, there are two main sections. The first section, 'Consensus mode', includes a 'Consensus type' dropdown menu set to 'Strict' and a 'Threshold' slider set to 100%. The second section, 'Export consensus', is highlighted with a red background and contains an 'Export to file' text box with the filename 'Aligned reads_consensus_3.txt', a 'File format' dropdown menu set to 'Plain text', a checkbox for 'Keep gaps' which is currently unchecked, and an 'Export' button at the bottom.

The following parameters are available:

Export to file - here you need to select path for the output file.

File format - format for the output file.

Keep gaps - set this parameter if the result consensus must keep gaps.

When you click on the *Export* button the consensus sequence will be exported into selected output file.