


# Selecting Sequence Region

## Selection in the Sequence View components


Holding the left mouse button down and moving the mouse, one can select a sequence region in the *Zoom View* or in the *Details View* components of the *Sequence View*. A sequence region that correspond to an annotation can also be selected by [double-clicking on the annotation](#) in the *Zoom View*, *Details View*, or *Annotations Editor*.

Selection of a sequence region is synchronized between different *Sequence View* components.


In the *Sequence Overview* the selection is displayed as a blue line:

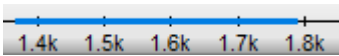

In the *Zoom View* a blue line with coordinates and the length of the selected region is shown:

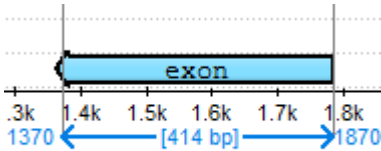


In the *Details View* a rectangle with a dashed border is shown around the selected region:

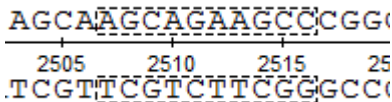


## Selecting a region using the mouse


It is possible to select a region of a sequence using the mouse in each component of the *Sequence View*.

In the *Zoom View* and the *Details View* the selection can be adjusted by dragging its left or right border:

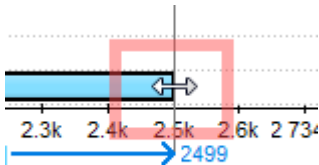

## Selecting a region using the dialog

If information about exact coordinates of a region is available, click the *Select sequence region* button  on the toolbar or *Select*→*Sequence region* in the context menu of the *Sequence View*. Input the coordinates in the *Region Selection* dialog:

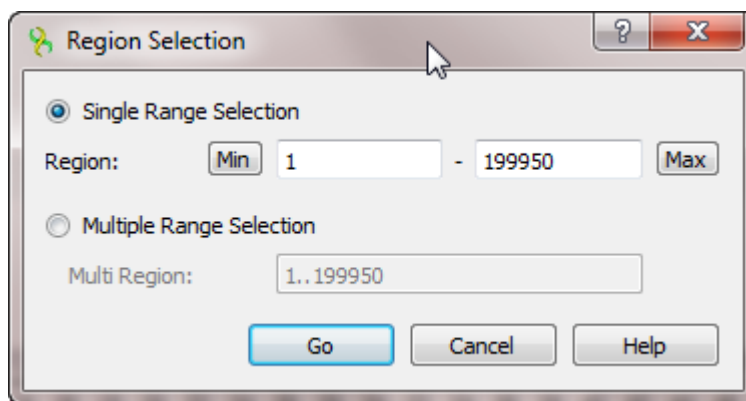


It is possible to input:

- *Single Range Selection* — a common region with a single range. Use the *Min* and *Max* buttons to automatically specify the beginning and the end of the sequence in the corresponding fields.
- *Multiple Range Selection* — a region that consists of several join elements, for example, exons of a gene. See "[The DDBJ/ENA/GenBank Feature Table Definition](#)" for details.

## Selecting a region around or between annotations

To select a region between two annotations, for example, an intron region between two exons:

- Click on the annotations, holding the *Ctrl/Cmd* key:

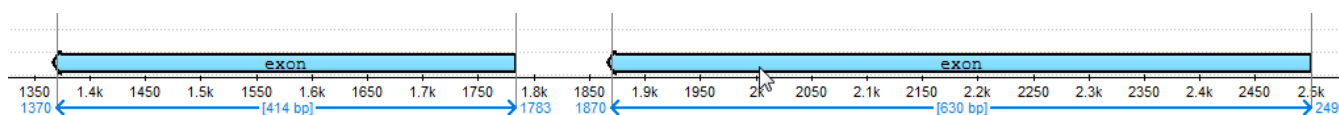

- Click *Select→Sequence between selected annotations* in the context menu of the *Sequence View*.

To select a region that contains the annotations in a single range, choose the *Sequence around selected annotations* item in the menu.