


Dashboard Window Components

Overview tab

"Workflow Task" widget

It contains:

- the workflow working time;
- the workflow running progress;
- the workflow task status: failed, finished successfully, running and etc.;

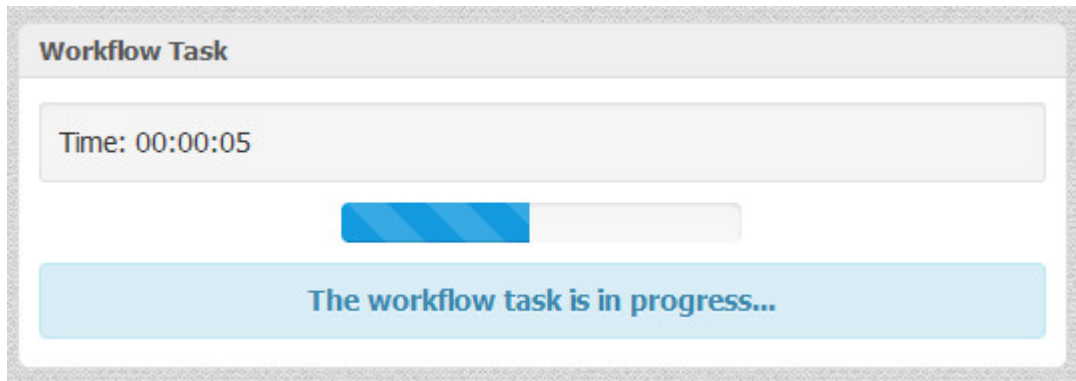

"Output Files" widget

It contains a table with the information about all created output files. The table columns are:

- clickable file name (with a help of the arrow on the right side of the file name you can open the file containing directory or open the file by operating system);
- the name of the workflow element that has produced the file;

Output Files	
File	Producer
muscle_alignment.aln	Write alignment
Open containing directory	

"Common Statistics" widget

It contains a table with common statistic for each workflow element in the workflow. The table columns are:

- name of the workflow element;
- time of the workflow element execution;
- the number of messages that has been retrieved;

Common Statistics		
Element	Elapsed time	Output messages
Align with MUSCLE	00:00:01	1
Read alignment	00:00:00	1
Write alignment	00:00:00	0

"Problems" widget

It contains a table with problems. The table columns are:

- problems type (warning, error and etc.)
- name of the element with problem
- error message

Problems		
Type	Element	Message
	Read Alignment	Unsupported document format

Input tab

"Parameters" widget

It contains a table with common statistic for each workflow element's parameter in the workflow. The table columns are:

- names of the workflow elements;
- names of the workflow parameters;
- values of the workflow parameters;
- clickable file name values of the workflow parameters (here you can open the file containing directory or open the file by operating system);

Parameters											
Read alignment											
Align with MUSCLE	<table><tr><th>Parameter</th><th>Value</th></tr><tr><td>Max iterations</td><td>-1</td></tr><tr><td>Mode</td><td>0</td></tr><tr><td>Region to align</td><td>Whole alignment</td></tr><tr><td>Stable order</td><td>True</td></tr></table>	Parameter	Value	Max iterations	-1	Mode	0	Region to align	Whole alignment	Stable order	True
Parameter	Value										
Max iterations	-1										
Mode	0										
Region to align	Whole alignment										
Stable order	True										
Write alignment											

External Tools tab

"External Tools" widget

It contains information about external tools. There are:

- names of the external tools;
- executable file of the external tool;
- arguments of the external tool;

The screenshot shows the 'External Tools' widget with two tool configurations:

- Find Peaks with MACS**
 - MACS run 1** (with a document icon)
 - Run info**
 - Executable file**: C:\Python27/python.exe
 - Arguments**
 - Error log**
- Create CEAS Report**
- Build Conservation Plot**
 - conservation_plot run 1** (with a document icon)
 - Run info**
 - Executable file**: C:\Python27/python.exe
 - Arguments**:

```
--phasdb=E:/UGENE/trunk/data/cistrome/phastCons/hg19
--height=1000
--width=1000 "-w 1000"
--title="Average Phastcons around the Center of Sites"
--bed-label=Conservation_at_peak_summits C:/Users/yalgaer/AppData/Local/Temp/ugene_tmp/p54244/ConservationPlot_tmp/13830
38905_0/Conservation_at_peak_summits.bed
```
 - Error log**

To copy external tool run string click the following button:

The screenshot shows a close-up of the 'External Tools' widget with the following configuration:

- Call Variants**
 - SAMtools run 1** (with a document icon highlighted by a red circle)
 - Run info**
 - Executable file**

With a help of the context menu of this widget you can copy selected text from the dashboard or copy all text of the active element:

Error log

INFO @ Tue, 29 Oct 2013 16:25:48:

ARGUMENTS LIST:

name = Default

format = AUTO

Copy selected text

Copy element content