


Exporting Sequence Image

Use a sequence toolbar *Export image* button to save a screenshot of the sequence:

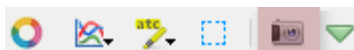

The Export Image dialog will appear where you should set name, location, export settings and format of the picture:

A screenshot of the 'Export Image' dialog box. The dialog has a title bar with a close button. It is divided into two main sections: 'Export settings' and 'Export to file'. Under 'Export settings', there are two columns. The 'Area' column has three radio buttons: 'Currently viewed', 'Zoomed annotations', and 'Sequence details' (which is selected). The 'Region' column has a dropdown menu showing 'Selected region' and two input boxes containing '31' and '75' separated by a hyphen. Under 'Export to file', there is a 'File name' field with a browse button (...), a 'Format' dropdown set to 'png', and a 'Quality' slider set to 80. At the bottom are three buttons: 'Help', 'Cancel', and 'Export'.

UGENE supports export to the BMP, JPEG, JPG, PNG, PPM, TIF, TIFF, XBM, XPM and SVG image formats. You can export currently viewed, zoomed annotations or sequence detailes areas. Also you can export whole sequence or custom region. Use the *Region* settings to do it.