


Copying Sequences

To copy current selection click the *Copy Copy selection* item in the *Actions* main menu or the context menu. The hotkey for this action is Ctrl-C.

To copy current selection in the different formats click the *Copy->Copy formatted* item in the *Actions* main menu or the context menu. The hotkey for this action is Ctrl+Shift+C. To select the copying format use the *General Options Panel* tab.

To copy one or several sequences do the following:

- Select the sequences in the *Sequence list* area;
- Select the *Copy Copy selection* context menu item in the *Sequence area* or use hot key combination. Note, that if you activate context menu in the *Sequence list* area you will lose your current selection.

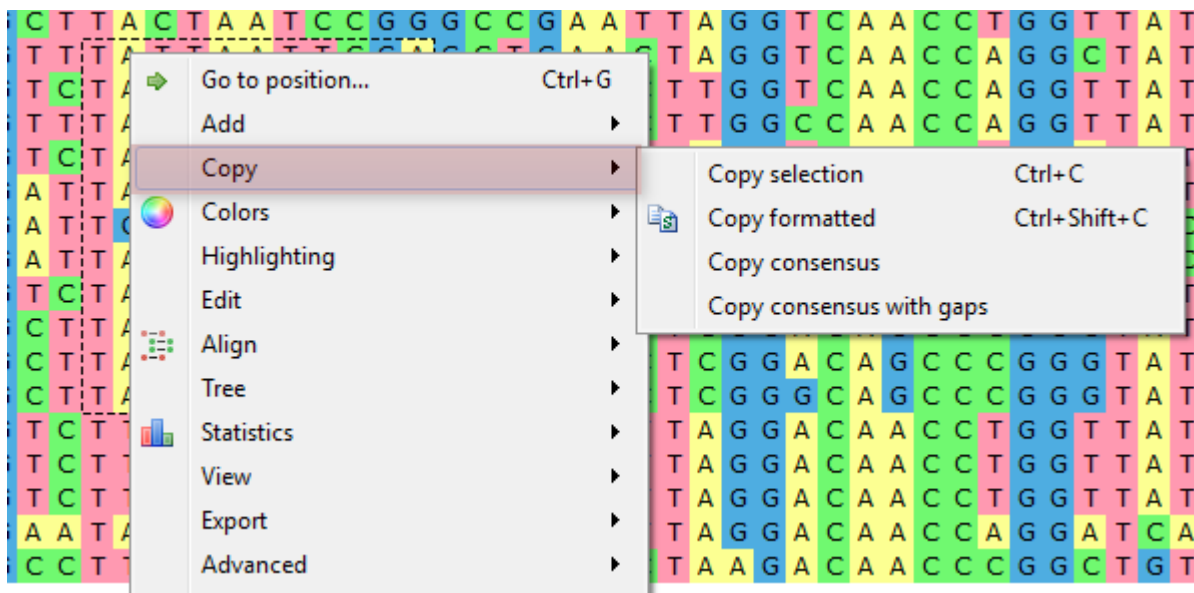

To copy consensus sequence use the *Copy Copy consensus* item.