
Creating Document

To create a new sequence file from a text, select the File New document from text main menu item.


The Create Sequence dialog appears:

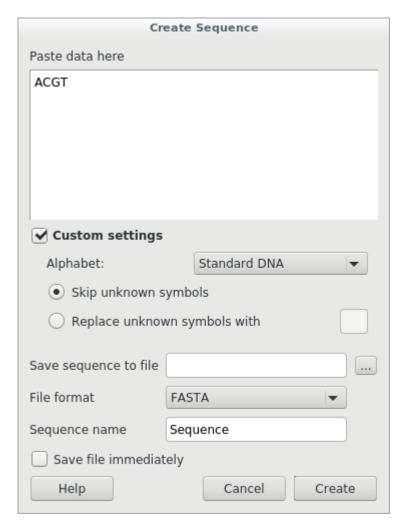

You can input the created sequence to the Paste data here field. You can type or paste sequences in FASTA or text format.

The following Custom settings are available:

Alphabet — here you can select the alphabet:

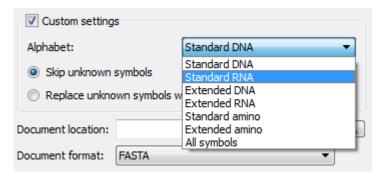

The following alphabets are available: Standard DNA, Standard RNA, Extended DNA, Extended RNA, Standard amino, Extended amino.

Skip unknown symbols / Replace unknown symbols with — you can select either to skip unknown input symbols or to replace them with the specified symbol.

Document location — a location of the created document.

Document format — a format of the created document. Currently available formats are FASTA and Genbank.

Sequence name — a name of the sequence in the created document.

Save file immediately — check this option if you want to save the document immediately after the *Create* button is pressed.

The created document will be added to the current project and opened in the Sequence View.