


Saving Graph Cutoffs as Annotations

To save graph cutoffs as annotations select the *Graph->Save cutoffs as annotations* item in the graph context menu. The following dialog will appear:

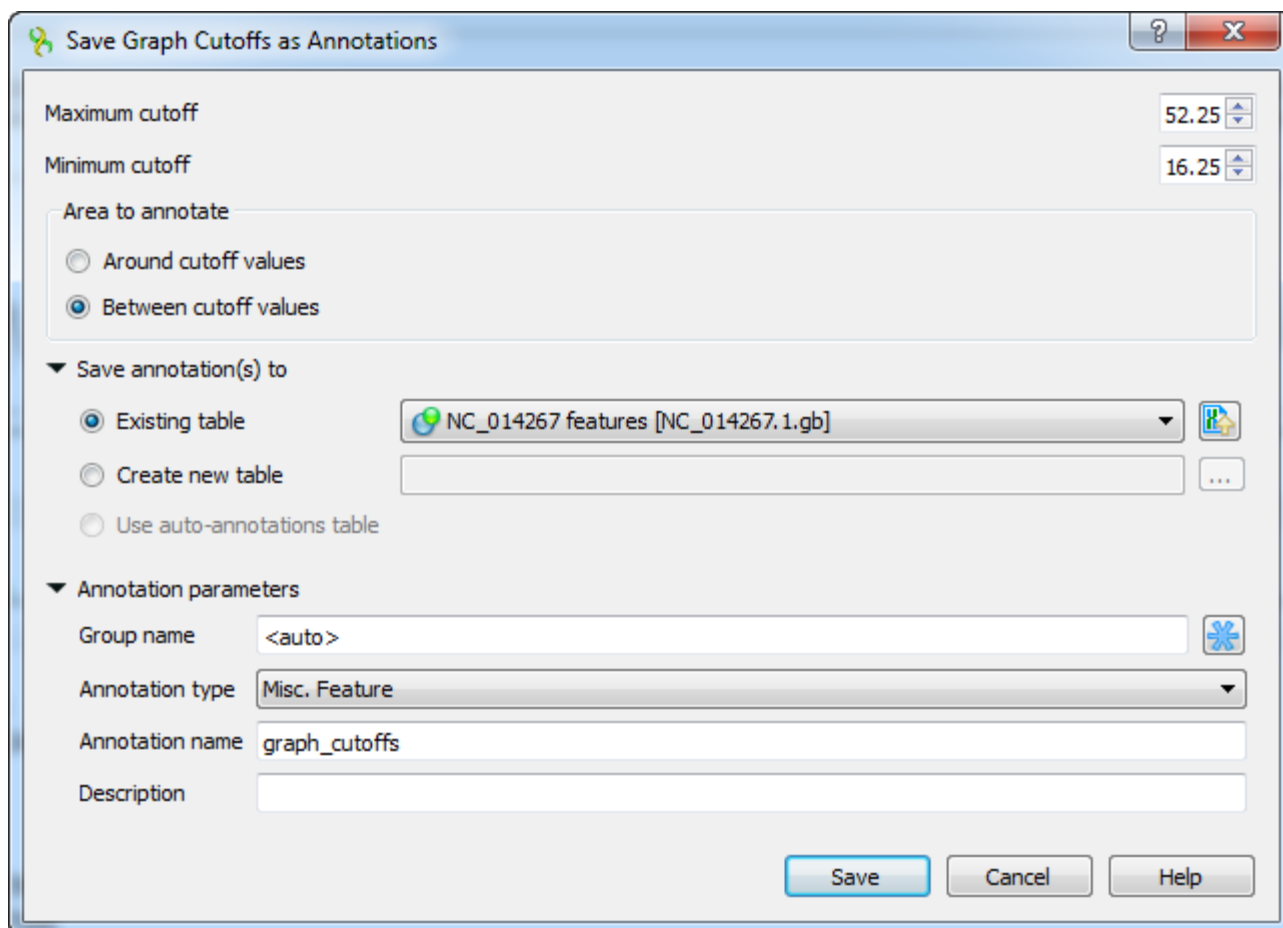

The dialog box is titled "Save Graph Cutoffs as Annotations". It contains the following fields and controls:

- Maximum cutoff:** A numeric input field with a value of 52.25 and a small up/down arrow icon.
- Minimum cutoff:** A numeric input field with a value of 16.25 and a small up/down arrow icon.
- Area to annotate:** A section with two radio buttons:
 - ☐ Around cutoff values
 - ☒ Between cutoff values
- Save annotation(s) to:** A section with three radio buttons:
 - ☒ Existing table: A dropdown menu showing "NC_014267 features [NC_014267.1.gb]" and a file icon button.
 - ☐ Create new table: An empty text input field and a "..." button.
 - ☐ Use auto-annotations table
- Annotation parameters:** A section with four fields:
 - Group name:** A text input field containing "<auto>" and a refresh icon button.
 - Annotation type:** A dropdown menu showing "Misc. Feature".
 - Annotation name:** A text input field containing "graph_cutoffs".
 - Description:** An empty text input field.

At the bottom right, there are three buttons: "Save", "Cancel", and "Help".

The following parameters are available:

Maximum cutoff - maximum cutoff value.

Minimum cutoff - minimum cutoff value.

Around cutoff values - saves the values around cutoffs values.

Between cutoff values - saves the values between cutoffs values.

In the *Save annotation(s) to* group you can set up a file to store annotations. It could be either an existing annotation table object, a new annotation table or auto-annotations table (if it is available).

In the *Annotation parameters* group you can specify the name of the group and the name of the annotation. If the group name is set to <auto> UGENE will use the group name as the name for the group. You can use the '/' characters in this field as a group name separator to create subgroups. If the annotation name is set to *by type* UGENE will use the annotation type from the *Annotation type*: table as the name for the annotation. Also you can add a description in the corresponding text field.

Select the parameters and click on the **Save** button. The corresponding annotations will be saved.