


Exporting Sequence from Alignment

To export one sequence from an alignment select the sequence in the sequence list or in the sequence area and use the *ExportExport selected sequence(s)* context menu item. The following dialog will appear:

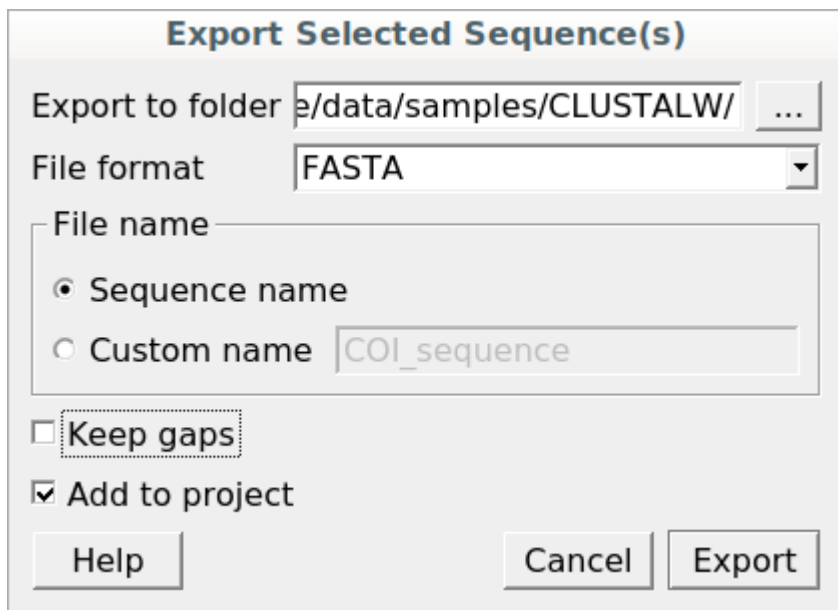

The dialog box is titled "Export Selected Sequence(s)". It contains the following fields and controls:

- Export to folder:** A text field containing the path "e:/data/samples/CLUSTALW/" followed by a browse button "...".
- File format:** A dropdown menu currently showing "FASTA".
- File name:** A section containing two radio buttons:
 - Sequence name:** The selected option.
 - Custom name:** A text field containing the text "COI_sequence".
- Keep gaps:** A checkbox that is currently unchecked.
- Add to project:** A checkbox that is currently checked.
- Buttons:** "Help", "Cancel", and "Export" buttons are located at the bottom of the dialog.

Here, it is possible to specify the result file location, file name, select a sequence file format, and define whether to keep or remove gaps ('—' chars) in the sequence and optionally add the created document to the current project.