

# PCR Primer Design for DNA Assembly

## Overview

The [polymerase chain reaction](#) (PCR) uses a pair of custom primers to direct DNA elongation toward each other at opposite ends of the sequence being amplified. These primers are typically between 18 and 24 bases in length and must code for only the specific upstream and downstream sites of the sequence being amplified. A primer that can bind to multiple regions along the DNA will amplify them all, eliminating the purpose of PCR.

A few criteria must be brought into consideration when designing a pair of PCR primers. Pairs of primers should have similar melting temperatures since annealing during PCR occurs for both strands simultaneously, and this shared melting temperature must not be either too much higher or lower than the reaction's [annealing temperature](#). A primer with a  $T_m$  (melting temperature) too much higher than the reaction's annealing temperature may mishybridize and extend at an incorrect location along the DNA sequence. A  $T_m$  significantly lower than the annealing temperature may fail to anneal and extend at all.

Additionally, primer sequences need to be chosen to uniquely select for a region of DNA, avoiding the possibility of hybridization to a similar sequence nearby.

Also, there are other parameters to be taken in account - [Gibbs free energy](#), self- and heterodimers availability, etc. All parameters are described below in details.

UGENE provides the "PCR Primer Design for DNA Assembly" feature only for nucleic sequences with the "Standard DNA" alphabet. To use it in UGENE open a DNA sequence and go to the "PCR Primer Design for DNA Assembly" tab of the Options Panel:

PCR Primer Design for DNA assembly

User primers

Forward user primer:

Reverse user primer:

Choose generated sequences as user primer's end

Filter:

|   | Sequence |
|---|----------|
| 1 | TTTGTGTC |
| 2 | TTTGTGGT |
| 3 | TTTGTGGG |
| 4 | TTTGTGGA |
| 5 | TTTGTGAG |
| 6 | TTTGTCTG |

Find reverse complement

Add sequence to a user primer's end:

Forward 5'

Reverse 3'

Parameters of priming sequences

Gibbs free energy: -40 kcal/mole  - 30 kcal/mole

Melting point: 57 °C  - 65 °C

Overlap length: 18 nt  - 25 nt

Parameters to exclude in whole primers

Gibbs free energy: -13 kcal/mole

Melting point: 20 °C

Max complementary motifs length: 3 nt

Select areas for priming search

Primer search areas for insert

Left area: 1  - 71

Right area: 355  - 426

Open the backbone sequence

Insert to backbone bearings

☒ 5' insert to 5' backbone ☐ 3' insert to 3' backbone

5' backbone length: 0  3' backbone length: 0

Other sequences in PCR reaction

| Fragment | Region |
|----------|--------|
|----------|--------|


There are the following parameters:

- User primers - the pair of forward and reverse primers user can check on the absence of unwanted connections (like hairpins, self- and hetero-dimers).
- Choose generated sequences as user primer's end - user primers could be extended with the 8-bases sequences on 5' and 3' ends. These sequences do not form too simple unwanted connections (for example, "AAAATTTT", obviously, has a hairpin, so there is no such a sequence in the list).
- Parameters of priming sequences - the range of melting temperature, Gibbs free energy and length primers SHOULD have.
- Parameters to exclude in whole primers - the minimum melting temperature, the maximum Gibbs free energy and the maximum base pairs length which may be present in the self- and hetero-dimers in the considered primer.
- Select areas for priming search - select areas to search forward and reverse primers in.
- Open the backbone sequence - choose a sequence to be added to 5' or 3' end of the result primer.
- Other sequences in PCR reaction - the list of sequence to check the unwanted connections with. The unwanted connections should fit to "Parameters to exclude in whole primers".

## Designing primers

Physical quantity

The designed primers should have the set Gibbs free energy, melting temperature and length. Also, these primers could have hairpins, self- and hetero-dimers, but the mentioned parameters of these found dimers should be limited. Look at the pictures:

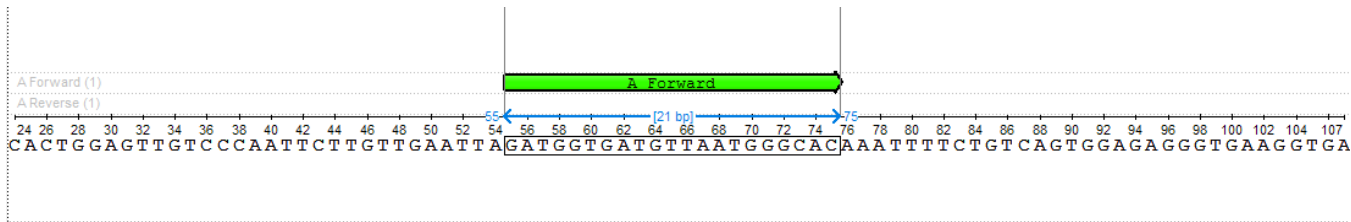

▼ Parameters of priming sequences

|                   |               |   |               |
|-------------------|---------------|---|---------------|
| Gibbs free energy | -40 kcal/mole | - | -30 kcal/mole |
| Melting point     | 51 °C         | - | 65 °C         |
| Overlap length    | 18 nt         | - | 25 nt         |


▼ Parameters to exclude in whole primers

|                                 |              |
|---------------------------------|--------------|
| Gibbs free energy               | -7 kcal/mole |
| Melting point                   | 20 °C        |
| Max complementary motifs length | 3 nt         |

The found primer "GATGGTGATGTTAATGGGCAC" has Gibbs free energy between -40 and -30 kcal/mol, melting temperature between 51 and 65 °C and length between 18 and 25 nucleotides. Also, there are no hairpins, self- and hetero-dimers in this primer, which has Gibbs free energy more than -7 kcal /mol, melting temperature less than 20 °C and base pairs length more than 3 nucleotides.

Regions

Forward primers are searched at the left region, reverse primers are searched at the right region. There is an amplified fragment between them:

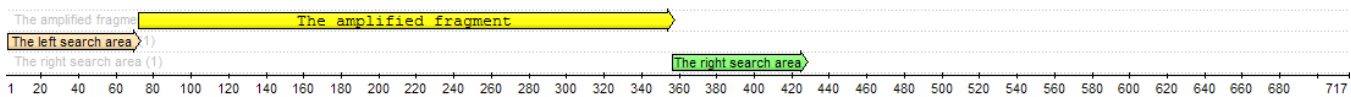

▼ Select areas for priming search


Primer search areas for insert

|             |     |   |     |                 |
|-------------|-----|---|-----|-----------------|
| Left area:  | 1   | - | 71  | Select manually |
| Right area: | 355 | - | 426 | Select manually |

The left search area is from 1 to 71 base, the right search area is from 355 to 426 base. The amplified fragment between them.

Backbone sequence

Backbone sequence is the sequence that is added to 5' end of the primer for the following assembly of the PCR product with the backbone molecule. This parameter is optional. You may set the backbone sequence, the sequence end you want to add backbone to and the length of the sequence from 5' and 3' ends to check if these sequences have unwanted hairpins, self- and hetero-dimers. On the pictures below you may see the differences between "5' insert to 5' backbone" and "5' insert to 3' backbone" ("insert" in this case is the same as "primer").

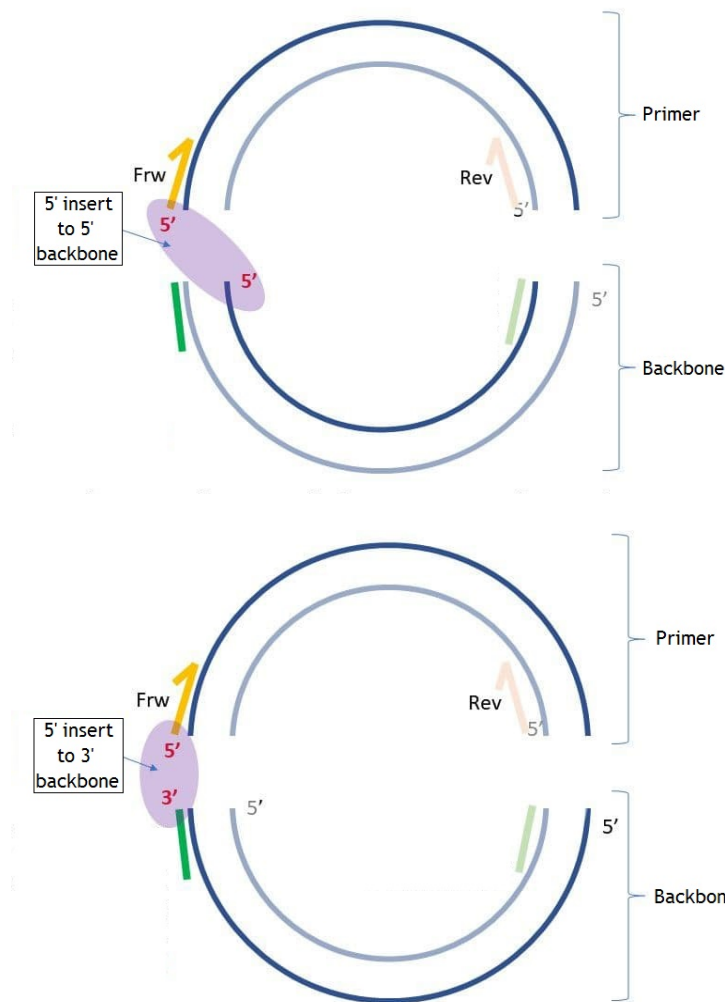

If you set "5' backbone length", then the number of bases you set from 5' end of the backbone sequence will be checked for unwanted connections. The same about "3' backbone length", but from the 3' end. If the checked end has hairpins, self- and homo-dimers you will see the following dialog:

- [Backbone details](#)

## Other sequences in PCR reactions

Sequences that are supposed to be in the reaction mixture. Its important primers do not have unwanted connections to these sequences, so it will be checked that there are no hairpins, self- or hetero-dimers between the resulting primer and any sequence from the "Other sequences in PCR reaction" file.

## Result

The search result is from 1 to 4 pairs of primers. The results should be as close as possible to the amplified fragment. That means, that the searching process starts from the edge of the amplified fragment to the 5' direction on the direct strand in case of the forward primer and to the 3' direction of the reverse-complementary strand in case of the reverse primer. The results are defined with special names:

- A - the pair of primers closest to the amplified fragment, corresponding to the entire set of parameters.
- B1 - the pair of primers closest to the amplified fragment, corresponding to the entire set of parameters if B2 and/or B3 exist.
- B2 - the pair of primers corresponding to the entire set of parameters, which has a 4-nucleotide intersection with B1.
- B3 - the pair of primers corresponding to the entire set of parameters, which starts immediately when B1 ends.

